




MEDPOR TITAN[®] Cranial Curve


Designed for cranial reconstruction

MEDPOR TITAN Cranial Curve

- The MEDPOR TITAN Cranial Curve has the strength of titanium with the flexibility and tissue in-growth¹ benefits of MEDPOR biomaterial
- Implants are packaged sterile and incorporate titanium mesh within a pre-contoured MEDPOR implant that is designed to mimic the native curvature of the cranium
- Two implant configurations are available:
 - Completely porous implant may allow for tissue in-growth on both sides
 - BARRIER option incorporates a non-porous sheet of polyethylene on the inferior aspect of the implant to prevent tissue in-growth from the dural side while allowing tissue in-growth from the porous, superior side
- Implants may be contoured intra-operatively with surgical scissors to fit the defect site² as needed
- The titanium mesh within the implant makes them easily identifiable on radiographs and CT scans
- May be fixated with any Stryker Universal Neuro self-drilling or self-tapping screws

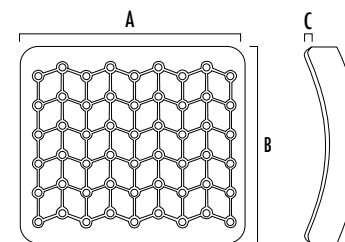

CAT#	DESCRIPTION	A	B	C
82019	TITAN Cranial Curve	38mm	48mm	0.85mm
82020	TITAN Cranial Curve - BARRIER	38mm	48mm	1.00mm

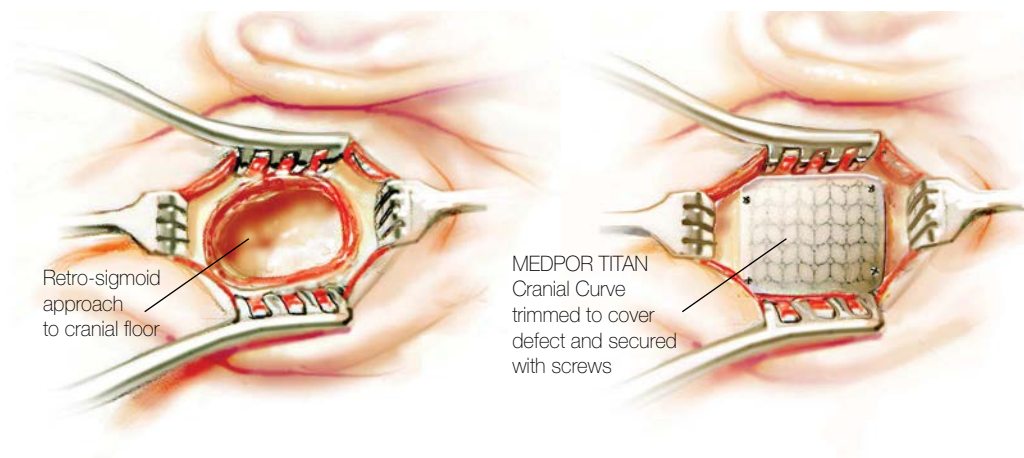

1. Liu JK, Gotfried ON, Cole CD, Dougherty,WR, CouldwellWT, "MEDPOR Porous Polyethylene Implant for Cranioplasty and Skull Base Reconstruction"Neurosurgery [April 2004]
2. Holck, D., Foster J., and Dahl T., "Custom Shaped Porous Polyethylene-Titanium Mesh Orbital Implants for Internal Orbital Floor/MedialWall Fracture Repair" ASOPRS 37th Annual Fall Scientific Syllabus, pp190, November 15-16, 2006

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the breadth of Stryker product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: **MEDPOR, MEDPOR TITAN, Stryker, Universal Neuro**. All other trademarks are trademarks of their respective owners or holders.

stryker[®]

Stryker Craniomaxillofacial
Kalamazoo,MI 49002 USA
t: 269 324 5346
toll free: 800 962 6558
f: 877 648 7114

www.stryker.com

Literature Number:
9410-400-256 Rev. None
UnDe/P.S.
Copyright © 2012 Stryker
Printed in USA